

Technical Data Sheet

Doc No – Ab/TDS/RB-43-II

Rev No : 02

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Test Tube Rack, PC

Valid for Cat. No. P20706, P20707, P20708, P20709, P20710, P20711, P20712, P20713, P20714, P20715

1. Description/Specification:

2. 1.1 Description:

Reusable Test Tube Rack PC, DNase, RNase and Pyrogen Free

1.2 Design:

Ergonomic design for easy handling. Holes are Keyed with molded-in- numbers and letters for easy identification.

1.3 Performance Specifications:

1.3.1 Capacity: P20706 (13 mm), P20707 (13 mm), P20708 (16 mm), P20709 (16 mm), P20710 (20 mm), P20711 (20 mm), P20712 (25 mm), P20713 (25 mm), P20714 (32 mm), P20715 (32 mm),

1.4 Material:

1.4.2 Complies with requirement of European Pharmacopoeia, 7th Edition (2011) & its supplements 7.5 (07/2012) monograph 3.2.2 Plastic container and closures for Pharmaceutical use, Confirm US FDA 21CFR.

1.4.3 Complies with the requirement of European Pharmacopoeia, 7th edition (2011) & its supplements 7.5 (07/2012) monographs 3.2.2 Plastic container and closures for Pharmaceutical use

1.4.4 The product has Passed the USP testing including class VI tests & has been assigned the FDA Drug master file DMF 12123

1.4.5 Material has Successful Passed the biological tests according to ISO 10993 external Communication devices for indirect blood contact for Prolonged period.

1.5 Color: White, Yellow, Blue



Shiv Ganga Industrial Estate,
Plot no. 1 Lakeshwari, Bhagwanpur,
Roorkee - 247667 Uttarakhand, India
An ISO 9001:2008 Company
A Division of Abdos

Regd. Office :
13A, Government Place (East)
Kolkata – 700069, India
CIN : U85100WB2008PTC123903

P : +91 (0) 1332 235277
E : info@abdoslaware.com
W : www.abdoslaware.com

Technical Data Sheet

1.6 Sterilization:

Sterile catalog number from this lot have been pre-sterilized and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation.

1.7 Quality control:

DNase, RNase, and Pyrogen Free

Testing of various characteristics in accordance with the valid Specification

1.8 Other information

Reusable, can be used in wide temperature range

1.9 Basic Features


1.9.1 Temperature range: -20⁰C to +121⁰C


1.9.2 Autoclavability: Yes, Autoclavable at 15 psi and 121⁰C for 20 minutes

1.9.3 Chemical resistance: See website.

Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

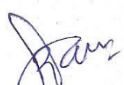
Quality Assurance Incharge
Abdos Labtech Pvt. Ltd.

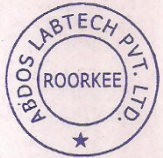


Shiv Ganga Industrial Estate,
Plot no. 1 Lakeshwari, Bhagwanpur,
Roorkee - 247667 Uttarakhand, India
An ISO 9001:2008 Company
A Division of Abdos

Regd. Office :
13A, Government Place (East)
Kolkata – 700069, India
CIN : U85100WB2008PTC123903

P : +91 (0) 1332 235277
E : info@abdoslabware.com
W : www.abdoslabware.com

Technical Data Sheet

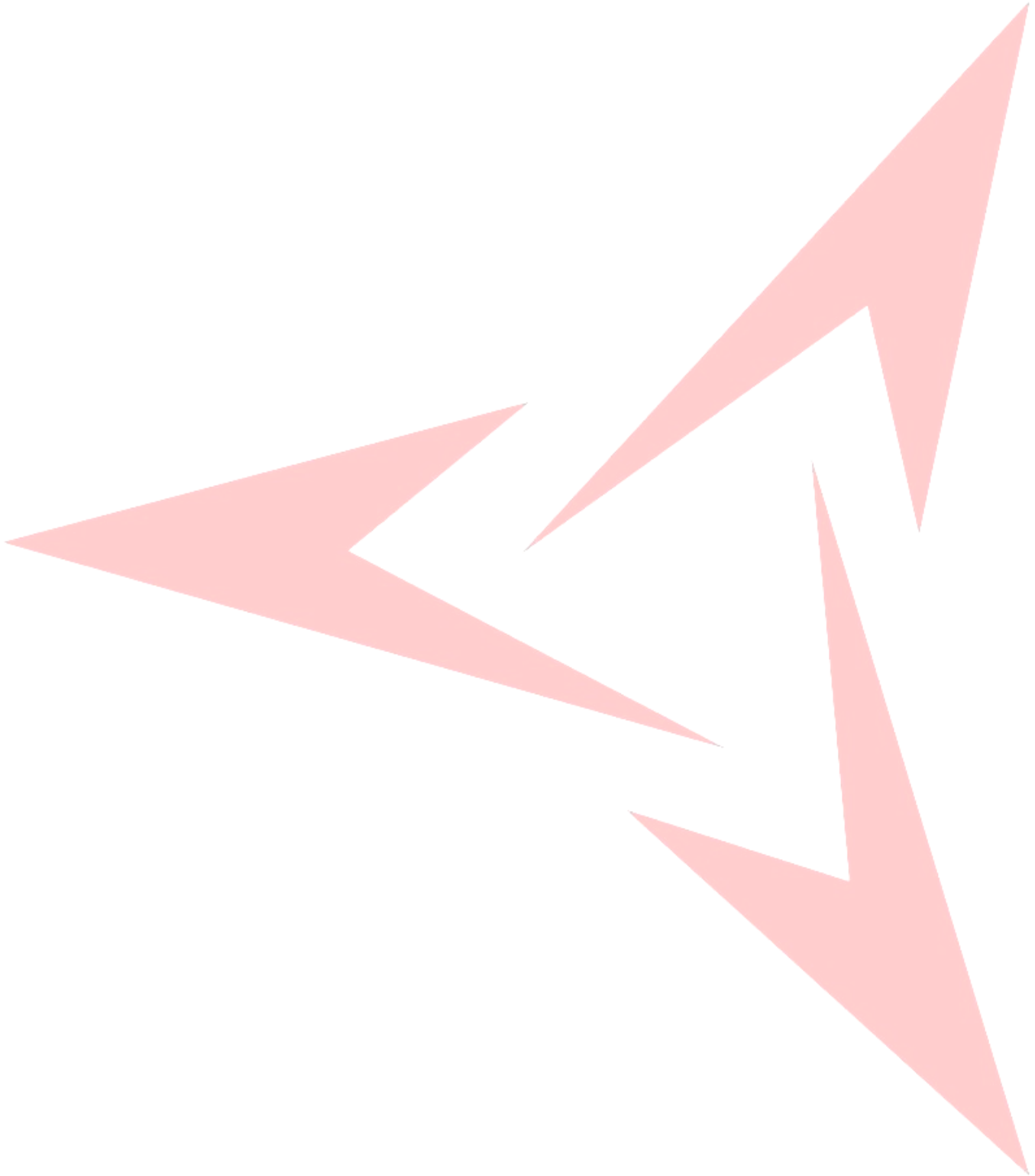

Shiv Ganga Industrial Estate,
Plot no. 1 Lakeshwari, Bhagwanpur,
Roorkee - 247667 Uttrakhand, India
An ISO 9001:2008 Company
A Division of Abdos

Regd. Office :
13A, Government Place (East)
Kolkata – 700069, India
CIN : U85100WB2008PTC123903

P : +91 (0) 1332 235277
E : info@abdoslabware.com
W : www.abdoslabware.com